

# 710001 Cell Scrapers

Specially designed Nest® cell scrape can be used to scrap and collect cells easily and effectively. Push the handle with a slight pressure to change the blade angle. Then push down the handle towards the bottom of the container. Rotate the handle slightly to make the blade twist to right direction.

## Features

1. Blade is made of HIPS and handle is made of PP
2. Free rotating blade, make the blade twist to right direction
3. Easy-tear sterile packing
4. Sterilized by E-beam, SAL=10<sup>-6</sup>
5. Non-pyrogenic and DNase/RNase free



## Description

Shelf Life : 3 Years after date of production ( ensure package is in good )

## Product Range

Cat.No.	Total Length ( mm )	Blade Width ( mm )	/Pack	/Case
710001	220	13	1	100

## Technical Drawing

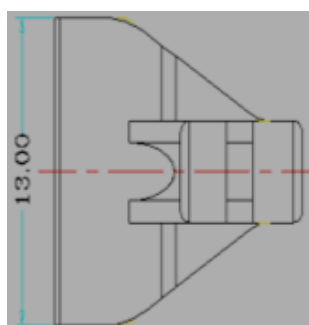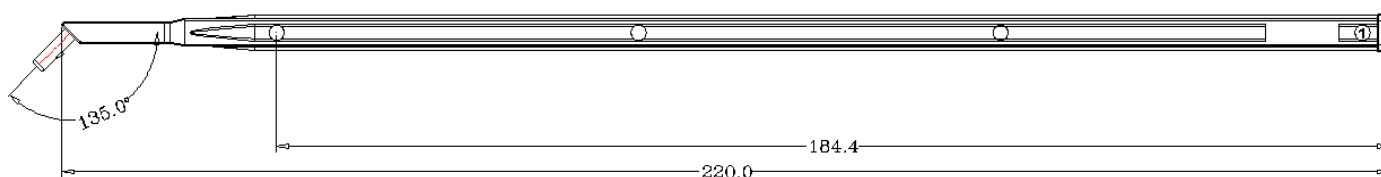

<b>NEST China</b>	Tel	: 0510-6878 8698 / 6878 8718	<b>NEST USA</b>	Tel	: +1-732-381-0268
	Email	: info@nest-wuxi.com		Email	: info@nestscientificusa.com
	Website	: www.cell-nest.com		Website	: www.nestscientificusa.com
	Add	: No.230, xida road, new district, Wuxi, Jiangsu, China		Add	: 1592 Hart St., Rahway, NJ07065, USA

For Research Use Only. Not intended for use in diagnostic or therapeutic procedures.

© Wuxi NEST Biotechnology Co., Ltd. Version 1, Jan. 10th, 2018