



# UVGI SANITIZING CABINETS

**Kills microorganisms and germs including viruses, bacteria and mold spores.**

**Our powerful germicidal UV-C sanitizes exposed surfaces of objects such as face masks, gloves, cell phones and tablets.**




## THE IDEAL UV SANITIZATION TOOL FOR HEALTHCARE PROTECTIVE WEAR




### UV SANITIZATION OF MASKS

Designed to make repetitive operations easy with **Auto-Repeat** functionality.

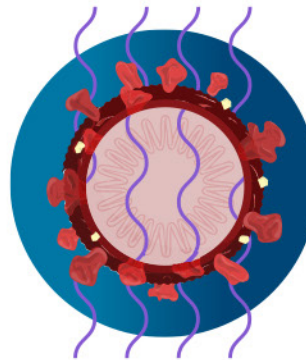

### HOW IT WORKS

UV-C light **disrupts** the genetic material inside **viruses, bacteria, and other microbes**, rendering them harmless.

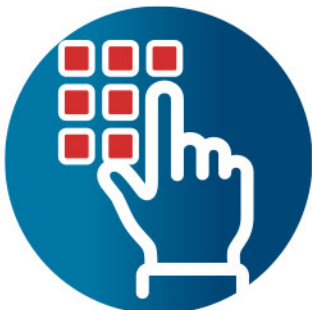

### EASILY PROGRAMMABLE

For sanitizing organism of choice: **continuous UV-C dosing** until full irradiation is performed with automatic shut-off.



### SAFE DURING OPERATION

Fully enclosed chamber **protects users** from accidental UV-C exposure due to opening of door.

# UV SANITIZATION DESIGNED FOR HEALTHCARE

## STANDARD

**HCL-1000**

Sanitize up to  
**3 N95 MASKS**

for full disinfecting without  
overcrowding



## LARGE

**HCL-1500**

Sanitize up to  
**6 N95 MASKS**

for full disinfecting without  
overcrowding



### Features:

- Embedded UV sensor inside chamber accurately measures UV-C output for adjusting dosage energy and time
- Programmable microprocessor circuitry automatically calculates irradiation time based on dosage energy entry, eliminating manual and guess work
- EMI/RFI protection
- Constructed of anodized aluminum for easy cleaning and durability

## PRODUCT SPECIFICATIONS:

<b>PRODUCT NUMBER:</b>	<b>HCL-1000</b> Standard Size	<b>HCL-1500</b> Large Size
<b>WAVELENGTH:</b>	254 nm	254 nm
<b>TUBES:</b>	(5) 8W SW - BLE-8T254	(6) 15W SW - BLE-1T155
<b>OVERALL HOUSING DIMENSIONS:</b> (W X H X D)	19.5 x 10.5 x 9 in (49.5 x 26.7 x 22.9 cm)	24 x 10.5 x 14 in (61 x 26.7 x 35.6 cm)
<b>EFFECTIVE INNER CHAMBER DIMENSIONS:</b> (W X H X D)	13.5 x 7 x 7.5 in (34.3 x 17.8 x 19.1 cm)	18.25 x 6.25 x 12.5 in (46.4 x 15.9 x 31.8 cm)
<b>DOOR DIMENSIONS:</b> (W X H)	12 x 6.25 in (30.5 x 15.9 cm)	16 x 6.25 in (40.6 x 15.9 cm)
<b>VIEWING WINDOW DIMENSIONS:</b> (W X H)	6 in x 2.5 in (15.2 x 6.4 cm)	6 x 2.5 in (15.2 x 6.4 cm)
<b>NET WEIGHT:</b>	17.5 lb (7.9 kg)	26 lb (11.8 kg)
<b>POWER:</b>	120V-60Hz-2A, 230V-50Hz-2A, 100V-50/60Hz-2A	120V-60Hz-3A, 230V-50Hz-3A, 100V-50/60Hz-3A