

500 ml Erlenmeyer Flask

Ideal for shaker culture applications.


Features

1. PETG Bottle & HDPE Cap
2. Sterilized by E-Beam
3. Nonpyrogenic
4. Packaged in sterile bag
5. Two types of cap style : vent caps & plug seal cap
6. Vent caps with 0.22um hydrophobic filters for gas exchange without contamination

Description

Cat.No. : 783001、783011
 Volume : 500ml
 Cap Style : Plug Seal Cap (783001)、Vent Cap (783011)
 Shelf Life : 3 Years after date of production (ensure package is in good)
 Sterile : Yes
 /Pack : 1
 /Case : 12

Technical Drawing

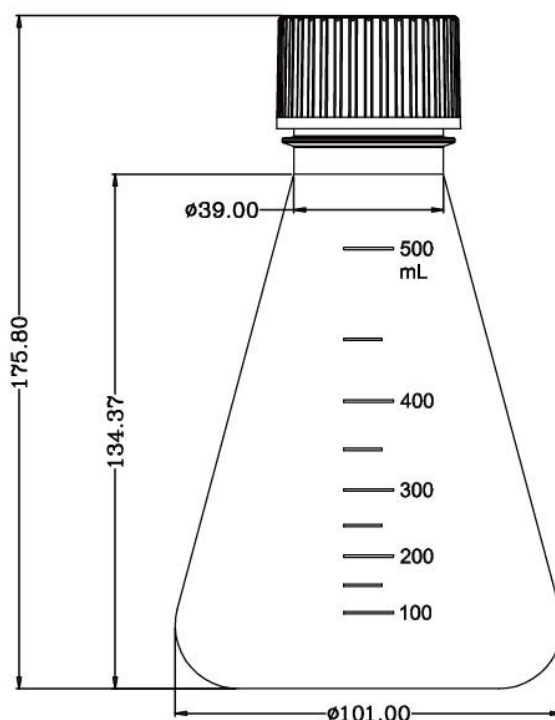

**NEST
China**

Tel : 0510-6878 8698 / 6878 8718
 Email : info@nest-wuxi.com
 Website : www.cell-nest.com
 Add : No.230, xida road, new district, Wuxi, Jiangsu, China

**NEST
USA**

Tel : +1-732-381-0268
 Email : info@nestscientificusa.com
 Website : www.nestscientificusa.com
 Add : 1592 Hart St., Rahway, NJ07065, USA

For Research Use Only. Not intended for use in diagnostic or therapeutic procedures.

© Wuxi NEST Biotechnology Co., Ltd. Version 1, Jan. 10th, 2018

Product Range

Cat.No.	Volume	Cap Style	Size (mm)			Sterile	/Pack	/Case
			Height	Bottleneck Diameter	Bottom Diameter			
781001	125ml	Plug Seal Cap	106.8	34	66	Yes	1	24
781011		Vent Cap						
782001	250ml	Plug Seal Cap	137.8	34	83	Yes	1	12
782011		Vent Cap						
783001	500ml	Plug Seal Cap	175.8	39	101	Yes	1	12
783011		Vent Cap						
784001	1000ml	Plug Seal Cap	213.3	39	127	Yes	1	6
784011		Vent Cap						

Recommendations

1. Recommended volume of culture medium : 30%-40% of bottle volume .
2. Suggested starting shaking speed of orbital Shaker : 75-125 rpm.

NEST China	Tel	: 0510-6878 8698 / 6878 8718	NEST USA	Tel	: +1-732-381-0268
	Email	: info@nest-wuxi.com		Email	: info@nestscientificusa.com
	Website	: www.cell-nest.com		Website	: www.nestscientificusa.com
	Add	: No.230, xida road, new district, Wuxi, Jiangsu, China		Add	: 1592 Hart St., Rahway, NJ07065, USA